



P.O. Box 219 • Batavia, Illinois 60510 • 1-800-452-1261 • flinn@flinnsci.com • Visit our website at: www.flinnsci.com

© 2002 Flinn Scientific, Inc. All Rights Reserved.

The Spinning Can

Lenz's Law Demonstration

Publication No. 10240

Introduction

Use Lenz's Law to "magically" spin an aluminum can without touching the can.

Science Concepts

- Lenz's Law
- Induced current

Materials

Aluminum can (with entire top removed)	Small bar magnet
Container (to float aluminum can)	Thread
Spill tray	Tap water

Safety Precautions

This demonstration does not present any unusual safety hazards. Care should be taken in removing the can covers.

Procedure

1. Locate an aluminum can (as heavy gauge as possible) with a diameter about 1" greater than the length of your bar magnet. An aluminum density catch bucket works well. (If the top of the can needs to be removed, do it carefully.)
2. Tie a thread tightly around the center of the bar magnet so that it can be spun in a circular fashion as shown in Figure 1.

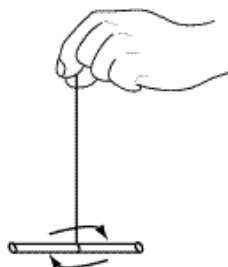


Figure 1. Practice spinning and balancing the magnet.

3. Locate a container (finger bowl, regular bowl, etc.) that is about 1" greater in diameter than the aluminum can.
4. Place the bowl in the center of a spill tray.
5. Fill the bowl to its very top with tap water until it can not hold any more water. A "bulging" bowl of water is ideal.
6. Carefully float the aluminum can open-side-up on top of the water. The final setup should look like Figure 2.
7. With the magnet hanging outside the can, start the magnet spinning on the hanging thread. Spin and balance the magnet so it is spinning in a level, even fashion. Then slowly lower the spinning magnet down inside the lip of the can, being careful to not touch the side of the can. (This might take a little practice and a steady hand.)

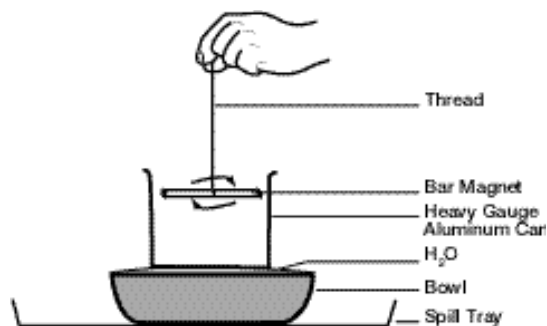


Figure 2

8. Watch what happens to the floating can. What direction does it spin compared to the magnet? What happens when the thread is completely unwound and the magnet spins in the opposite direction?

Discussion

In 1834 Henrich Lenz put forward his Law of Induction: *An induced current in a closed conducting loop will appear in such a direction that it opposes the motion that produced it.*

This "induced opposition" (often called Eddy currents) has been demonstrated in many ways since Lenz first formulated the principle. A classic demonstration has been to drop a strong magnet through a conducting tube to show how the created Eddy currents "oppose" the pull of gravity. The induced opposition causes the magnet to fall *very slowly* through the tube. It is a great demonstration to show the relationship between Lenz's Law and the Law of Conservation of Energy.

The spinning aluminum can in this demonstration illustrates the same principles. The can will spin in this demonstration in response to the induced current created between the spinning magnet and the aluminum can. The can is actually being repelled in response to the induced current from the spinning magnet. The "phase shift" is very dramatic when the induced field direction is reversed and the magnet starts spinning in the opposite direction. The can "magically" reverses its spinning direction too.

Tips

- Practice spinning and balancing the magnet prior to the demonstration as shown in Figure 1.
- The aluminum can should be as thick walled as possible. A soda can does **not** work well and has sharp edges when the entire top is removed. An aluminum catch bucket works well.
- This demonstration can be very dramatic if performed on an air table instead of on floating water.
- The right combination of sizes of parts is what makes this demonstration successful. The can needs to float on a high meniscus of water and therefore the bowl should be only slightly bigger in diameter than the can. If floating properly, the can should not touch the sides of the bowl. (A tea cup of the right size works well.)

Disposal

Pour water down the drain and save all of the other materials.

Reference

Dindorf, W. Lenz's Law in the Kitchen. *The Physics Teacher*. Vol. 37, May, 1999.

Lenz's Law products are available from Flinn Scientific, Inc.

Catalog No.	Description	Price/Each
AP4698	Eddy Current Experiment Kit	Consult Your Current <i>Flinn</i> Catalog/Reference Manual.
AP9261	Bar Magnet, Cobalt Steel	
AP4626	Catch Bucket, Aluminum	